

# It's a New Era in Power Mobility...



The Jazzy Air® is a new era in power mobility. It's social mobility. Whether at home, or on the town, the unique appeal of Jazzy Air will elevate your independence and social confidence in style.

- **Elevates** in just **16 seconds** to the primary standing height position allowing quicker access to your environment
- **Safely** drives **3.5mph** while elevated, enabling you to socialize with others at walking speed
- **Patented Active-Trac®** suspension provides **maximum stability** in both elevated and non-elevated positions

RECIPIENT OF:



## It's **Social** Mobility.

**Pride®**  
Mobility Products Corp.

# Jazzy air®

# Technical Specifications

Weight Capacity	300 lbs.
Maximum Speed <sup>1</sup>	Up to 4.0 mph (non-elevated position) Up to 3.5 mph (elevated position)
Ground Clearance <sup>2</sup>	2.42" at footplate 1.00" at motor bracket
High-Back Seat to Floor Range	23" - 32.9"
Turning Radius <sup>2</sup>	22.5"
Overall Size <sup>2,6</sup>	Length: 42.23" with foot rigging; 38.6" w/out foot rigging Width: 23.4" at the wheels 25.5" with armrests
Drive Controller	PG VR2
Drive Wheels	10" flat free
Caster Wheels	6" flat free
Suspension	Active-Trac <sup>®</sup> w/ Mid-Wheel 6 <sup>®</sup> Technology
Range <sup>1,4</sup>	Up to 15.29 miles w/ U1 batteries Up to 18.4 miles w/ 40 Ah batteries
Component Weights <sup>5</sup>	Base: 225 lbs. w/ U1 batteries 245 lbs. w/ 40 Ah batteries Batteries: U1-24.5 lbs. 40 Ah-32.2 lbs.
Drivetrain	Two motor, Mid-Wheel-6 drive
Brakes	"Intelligent Braking" (electronic regenerative, disc park brakes)
Batteries <sup>3</sup>	(2) 12 V, U1 (2) 40 Ah
Battery Charger	3.5 amp, off-board
Max. Allowable Incline	6°
Max. Angle Capacity	6°
Maximum Obstacle Climbing Ability	Lowered Position: 2.2" Elevated Position: 1.8"
Warranty	5 year limited (frame) 13 month limited (electronics) 13 month (drive motors)



## Features


- Achieve 10" of power adjustable seat height in 16 seconds
- Elevate or lower seat height while on the move with the "Air" button or PG VR2 joystick drive controller
- Seat and integrated footplate design provides maximum comfort in both elevated and non-elevated positions
- Angle adjustable footplate with three height adjustments
- Choice of 16" x 16", 18" x 18" and 20" x 20" contoured high-back reclining sport seat
- Angle, height & depth adjustable armrest
- Bright LED marker lights
- 60" lap belt

## Accessories

- Rear basket
- Weather cover
- Cup holder

## Shroud Colors

NOW Available:

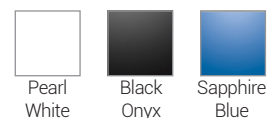

Coming Soon:



## Compatible Lifts

- Backpacker<sup>®</sup> Plus
- Outlander Full Platform
- Backpacker<sup>®</sup> AVP
- Outlander LP

## Options

- Lap belt: 50", 70", 80", 90"
- PG GC programmer

<sup>1</sup> Varies with base model, user weight, terrain type, battery amp hour rating (Ah), battery charge, battery condition, motors, controller type, tire type, and tire condition. This specification is subject to a variance of +10%, -5%.

<sup>2</sup> Due to manufacturing tolerances and continued product improvement, this specification can be subject to a variance of (+ or -) 3%.

<sup>3</sup> AGM or Gel-Cell type required. See VI. "Batteries and Charging."

<sup>4</sup> Tested in accordance with ANSI/RESNA, WC Vol. 2, Section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

<sup>5</sup> Battery weight may vary based on manufacturer.

<sup>6</sup> Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.

NOTE: This product conforms to all applicable ANSI/RESNA, ISO 7176 Series, and EN12184 standards. All specifications subject to change without notice.

AUTHORIZED PRIDE<sup>®</sup> PROVIDER

