

## MODEL

Oasis Collection	LC-580iM	LC-580iL
Position	True Infinite Position Zero Gravity Trendelenburg Chaise Lounger 5'4" - 6'	True Infinite Position Zero Gravity Trendelenburg Chaise Lounger 5'9" - 6'2"
Weight Capacity	375 lbs	375 lbs
Seat-to-Floor	20"	21"
Seat Depth	21"	21"
Seat Width	20"	21"
Top of Back to Seat	26"	29"
Back Style	Blown Fiber Biscuit Back	Blown Fiber Biscuit Back
Fully Padded Chaise	Yes	Yes
Seating Type	Pocketed Coil Springs w/Memory Foam Topper	Pocketed Coil Springs w/Memory Foam Topper
Heat & Massage	Yes	Yes
Distance From Wall	20"	21"
Overall Width	39"	39"
Chair Weight	130 lbs.	132 lbs.
Warranty	7-Yr. Pro-rated*	7-Yr. Pro-rated*

- Batteries not included.
- Crypton® Fabric offers a five-year warranty on the moisture barrier. Does not cover punctures in the fabric.

NOTE: Measurements may vary +/- 1". Specifications subject to change without prior notice.

Recommended HCPCS: E0627—Seat lift mechanism, electric, any type DME providers are responsible for determining the appropriate billing codes when submitting for insurance reimbursement. HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. Payer coding, coverage and bundling guidelines may apply.



## FABRIC SELECTION

**QS** **SARATOGA**

CASHMERE CHARCOAL EGGPLANT GODIVA NAVY

**LEXIS STA-KLEEN®** sta+kleen. PERFORMANCE FABRIC

BLACK MUSHROOM BURGUNDY CHESTNUT

**ULTRALEATHER™** Uf

BUFF FUDGE GARNET PECAN CHARCOAL

**ARIA<sup>2</sup>** CRYPTON. | super fabric

ESPRESSO COOL GREY RED LAZULI SAND

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice



Look for the Quick Ship Indicator for the shortest lead times available.

**Pride**  
Mobility Products Corp.

(US) 800-800-8586 • (Canada) 888-570-1113

pridemobility.com |

OASIS\_MEMORY\_CONTROL\_040419

PRELIMINARY

**Oasis**  
COLLECTION  
A New Level  
of Comfort




**POWER RECLINERS**

AMERICA'S #1 POWER LIFT RECLINER

pridemobility.com |

PRELIMINARY

*With the memory button hand control, your*  
**PERFECT POSITION**  
*is always at your fingertips!*

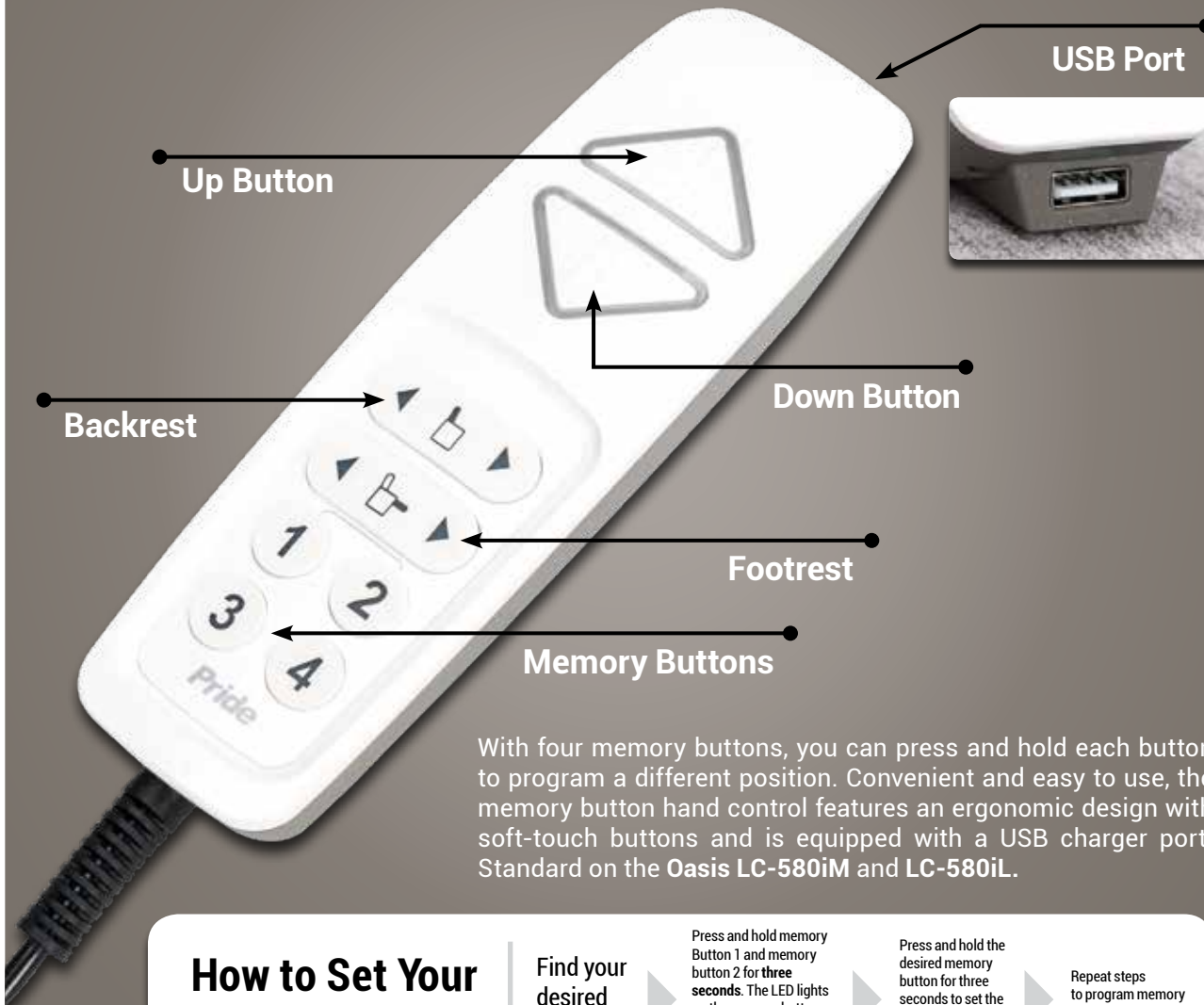

With four memory buttons, you can press and hold each button to program a different position. Convenient and easy to use, the memory button hand control features an ergonomic design with soft-touch buttons and is equipped with a USB charger port. Standard on the Oasis LC-580iM and LC-580iL.

### How to Set Your Memory Buttons

Find your desired position

Press and hold memory Button 1 and memory button 2 for **three seconds**. The LED lights on the memory buttons will flash to indicate that the hand control is in "set" mode.

Press and hold the desired memory button for three seconds to set the memory position. The LED lights will stop flashing.

Repeat steps to program memory buttons 2-4.



Shown in Cashmere

The Oasis brings a new level of comfort to Pride® Power Lift Recliners. Overstuffed padding provides elegant comfort and design, while contoured memory foam seating and pocketed coils deliver individualized support and relaxation. The Oasis provides true-infinite positioning and a standard memory button hand control.

- Premium furniture grade laminate/hardwood frames for added strength and durability
- Sinuous no-sag back springs combine durability with just the right amount of flexibility for personalized comfort
- Specially designed arms afford a luxurious feel and firm support
- Industry-exclusive quad pocket design offers extra storage – easy access top pockets provide convenient storage for your hand control and other devices
- Removable filling allows for custom lumbar support
- Standard lumbar pillow for additional support while relaxing

#### Options

- Deluxe heat and massage
- Residential cleaning kit
- Head pillow

**Oasis**  
COLLECTION