



SPECIFICATIONS

Overall Length (unfolded)	39 Inches
Overall Length (folded)	26 Inches
Overall Width	18 Inches
Ground Clearance	2.5 Inches
Maximum Load	320 Pounds
Range*	14.3 Miles
Total Weight (excluding battery)	51.7 Pounds
Turning Radius	35.5 Inches
Maximum Speed	4 MPH
Li-ion Battery Strength	10.5 Amp

*Range may vary with weight, terrain temperature, battery condition, etc. and is subject to manufacturing tolerances. Charge your battery when the indicator is red instead of letting the power run out for longer life of your battery.


LUGGIE ELITE FEATURES

- Alloy aluminum chassis frame
- Folds compactly for ease of transport & storage
- Super light weight, lifts easily into trunk of any car
- Battery removes easily
- Adjustable steering height and angle
- Comfortable seat with two height adjustments
- State of the art, speed control for safe operation
- Rear wheel drive train via direct trans axle
- 24 Volt DC permanent magnet motor
- Intelligent braking system (electronic, regenerative disc brakes)
- Plug-in battery charging system
- Free wheel facility



LED Battery Indicator is located on the top console of the Luggie scooter. As the battery is used, the LED indicator will move from green (fully charged) to red (low power).

Handlebars are extendable by one inch on either side for added comfort, flexibility and maneuverability.

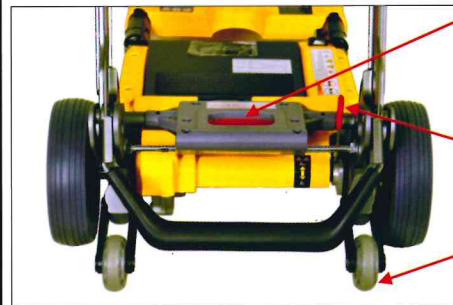

Seat Lock Lever located in the rear allows for 2 different height adjustments for driving comfort & suitability.

Free Wheel Lever puts the Luggie in free wheel mode allowing you to push with ease.

Anti-tipping Wheels is located under the rear bumper.



Sling Bag can be hung on the seat back as pictured or wore on the body. It comes with an orange reflector for driver safety.

Available Colors for the Luggie Elite are pictured. Colors shown left to right: blue, champagne, red, black, yellow and pink.

freerider USA

Luggie freerider