



## SPECIFICATIONS

Overall Length (unfolded)	39 Inches
Overall Length (folded)	26 Inches
Overall Width	18 Inches
Ground Clearance	2.5 Inches
Maximum Load	250 Pounds
Range*	11 Miles
Total Weight (excluding battery)	51.7 Pounds
Turning Radius	35.5 Inches
Maximum Speed	4 MPH
Li-ion Battery Strength	8.5 Amp


\*Range may vary with weight, terrain temperature, battery condition, etc. and is subject to manufacturing tolerances. Charge your battery when the indicator is in the red zone instead of letting the power run out, for longer life of your battery.

## LUGGIE FEATURES

- Alloy aluminum chassis frame
- Folds compactly for ease of transport & storage
- Super light weight, lifts easily into trunk of any car
- Battery removes easily
- Adjustable steering height and angle
- Comfortable seat with two height adjustments
- State of the art, speed control for safe operation
- Rear wheel drive train via direct transaxle
- 24 Volt DC permanent magnet motor
- Intelligent braking system (electronic, regenerative disc brakes)
- Plug-in battery charging system
- Magnetic brake release



**Battery Indicator** is located on the top console of the Luggie scooter. As the battery is used, the battery meter indicator needle will move from green (fully charged when needle is at far right) to red (low power).

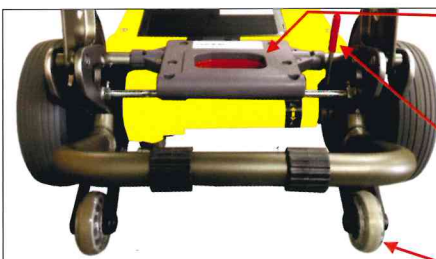

**Seat Lock Handle** located in the rear allows for 2 different height adjustments for driving comfort & suitability.

**Magnetic Brake Release** puts the Luggie in free wheel mode allowing you to push with ease.



**Anti-tipping Wheels** is located under the rear bumper.

**Sling Bag** can be hung on the seat back as pictured or wore on the body. It comes with an orange reflector for driver safety.

### Optional:

- Armrest
- Battery Bag
- Extra Battery
- Home Docking Station
- Sling Bag
- Soft Storage Case
- Suitcase

For the **Available Colors** for the Luggie contact your authorized Luggie reseller nearest you.

**FreeRider<sup>®</sup>USA**  
Luggie FreeRider