

VivaLIFT!®

POWER RECLINER COLLECTION¹

THE MOST
COMFORTABLE
RECLINER
IN THE WORLD!

THE 4-MOTOR RECLINER

NOW
AVAILABLE

TRUE-INFINITE
POSITIONING AND
MEMORY TOGGLE
HAND CONTROL!



Shown in Astro Brown Fabric



COUNTLESS
POSITION
COMBINATIONS

TRANQUIL 2¹

The **Tranquil 2** combines quiet comfort with gorgeous fabrics, for a rich and upscale look. Enhanced recline positioning to raise your feet above your heart.

The Tranquil 2 has been upgraded with an ergonomically designed memory toggle remote that features two programmable memory position buttons along with up/down/recline functions for customized comfort levels.



Scan here for
more information
on **VivaLift!**®









PRIDE®
MOBILITY
LIVE YOUR BEST®

TRANQUIL 2¹

VivaLIFT!
POWER RECLINER COLLECTION¹

COLLECTION	TRANQUIL 2 ¹	TRANQUIL 2 ¹	TRANQUIL 2 ¹	TRANQUIL 2 ¹
MODEL	PLR-935PW	PLR-935S	PLR-935M	PLR-935LT
POSITION	True-Infinite	True-Infinite	True-Infinite	True-Infinite
HEIGHT RANGE	5'3" and below	5'3" and below	5'4"-5'9"	5'10" and above
SIZE	Petite Wide	Small	Medium	Large/Tall
WEIGHT CAPACITY	400 lbs.	400 lbs.	400 lbs.	400 lbs.
SEAT-TO-FLOOR	18"	18"	19.5"	21.5"
SEAT DEPTH	18.5"	18.5"	20.5"	22.5"
SEAT WIDTH	24"	20"	22"	23.5"
TOP OF BACK TO SEAT	24.5"	25"	27"	30"
FULLY PADDED CHAISE	Yes	Yes	Yes	Yes
FOOTREST EXTENSION	Standard	Standard	Standard	Standard
SEATING TYPE	Total Comfort	Total Comfort	Total Comfort	Total Comfort
CHAIR HEIGHT	38"	40"	42"	47"
DISTANCE FROM WALL	12.5"	13"	16"	16"
OVERALL WIDTH	36"	33.75"	37.5"	38"
CHAIR WEIGHT	143 lbs.	137 lbs.	145 lbs.	155 lbs.
WARRANTY	7-Yr. Pro-rated	7-Yr. Pro-rated	7-Yr. Pro-rated	7-Yr. Pro-rated

INDEPENDENT FUNCTIONS

-  1. The power headrest lifts your head, neck and shoulders for maximum support.
-  2. The upward lift offers convenient stand assist for getting in and out of the recliner.
-  3. The power lumbar enables you to personalize your comfort.
-  4. The footrest extension allows you to relax and stretch out completely.
-  5. The backrest moves up and down independently to help you find the perfect position.
-  6. **NEW!** Enhanced true-infinite position recline enables the backrest and footrest to move independently so that a variety of positions can be reached including raising your feet above your heart.

STANDARD FEATURES

- 7. NEW!** Ergonomically designed toggle remote has two programmable memory position buttons. Performs the up/down, recline functions. Additionally, it has a convenient USB charger port that allows you to charge your smart device.
- 8.** Lithium battery backup ensures your recliner is operational when the power goes out.

1. VivaLift!® Power Recliners are FDA Class II Medical Devices designed to aid individuals with mobility impairments.


Note: Measurements may vary +/- 1". Specifications subject to change without prior notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice.



ES = Indicates shortest available lead time

Recommended HCPCS: E0627-Seat lift mechanism, electric, any type DME providers are responsible for determining the appropriate billing codes when submitting for insurance reimbursement. HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. Payer coding, coverage and bundling guidelines may apply.

STANDARD FABRIC: **ASTRO**

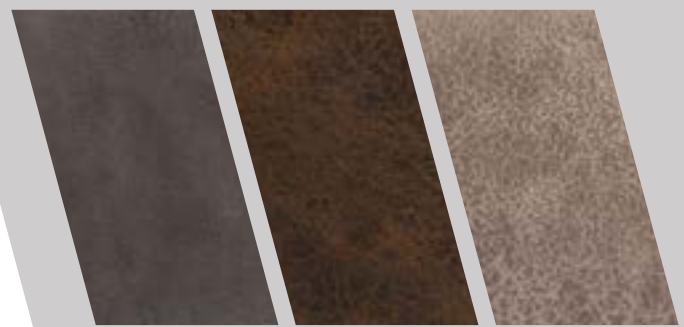

GREY

BROWN

MUSHROOM

(US) 800-800-8586 / (CAN) 888-570-1113 / pridemobility.com

PRIDE
MOBILITY
LIVE YOUR BEST®