

Model	440 HD PLUS Manual Folding
Overall Dimensions (L x W x H)	38" x 26" x 36"
Folded Dimensions (L x W x H)	29" x 16" x 26"
Seat Height (from ground)	20"
Armrest Height (from seat)	8"
Backrest Height (from seat)	18" - 24"
Max Loading Weight	440 lbs
Net Weight without Battery	61.5 lbs
Net Weight with 12AH Battery	68.5 lbs
Battery	24V 12AH or 24V 26AH
Battery Charge	AC 100-240, 50hz, 1.5-5A
Battery Charge Time	6 to 8 hours
Range with 12AH Battery	Up to 12 miles
Range with 26AH Battery	Up to 25 miles
Airline and Cruise Approval	Only with 24V 12AH Battery
Max Speed	4 MPH
Motor	2 * 250W (Total 500W) Brushless
Controller	24V 40A
Climbing Angle Max Degrees	9°
Max Safe Slope & Static Stability	6°
Brake	Intelligent electromagnetic brake
Turning Radius	≤47"
Reversing Width	≤49"
Front Wheel Diameter	8"
Rear Wheel Diameter	12"
Seating Area Width	20.5"
Cushion Size (width x depth)	19.7" x 18"
Controller Type	Detachable All-Direction Joystick
Warranty	Frame: 1 Year Drive Train: 1 Year Electronics: 1 Year Batteries: 1 Year

COMFYGO
mobility

makes life easier

comfygomobility.com

United States of America & Canada

- info@ComfyGoMobility.com
- +1 844 755 5055
- Visit comfygomobility.com/pages/contact or scan the QR code with your smartphone to fill out the help form.



Europe

- info@ComfyGoMobility.eu
- +31 (0) 20 809 7031
- Visit comfygomobility.eu/pages/contact or scan the QR code with your smartphone to fill out the help form.

COMFYGO
mobility



Loads
440lbs
Heavy Duty

makes life easier




MAJESTIC 440 HD PLUS
Remote Controlled Heavy
Duty Electric Wheelchair

FITS ON THE LEFT OR RIGHT ARMREST



360° Joystick

The wheelchair is equipped with a 360-degree universal intelligent joystick that offers easy control and maneuverability. The joystick includes a power on/off button with a power indicator light, speed selection buttons, and a horn for added convenience and safety.

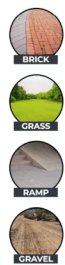

9° Climbing Gradient

intelligent electromagnetic brake



Customize Your Wheelchair

MAJESTIC 440 HD PLUS wheelchair comes with a raisable backrest and two battery sizes to choose from. You can also choose from four frame colors and four cushion colors to create your own unique style.

Airline & Cruise Approved

The MAJESTIC 440 HD PLUS electric wheelchair comes with storage bag under the seat and two battery size options. The 440 HD PLUS comes with a standard 12AH lithium-ion battery that is airline and cruise approved and offers up to 12 miles of driving distance. The optional 26AH lithium-ion, long-range battery is not airline or cruise approved but can travel up to 25 miles. Each battery can be charged while in the wheelchair or charged independently.



MAJESTIC 440 HD PLUS Remote Controlled Heavy Duty Electric Wheelchair

The MAJESTIC 440 HD PLUS wheelchair is a beautifully designed folding electric wheelchair for travel. This transport motorized wheelchair is lightweight, airline friendly. This is a sturdy wheelchair and features a weight capacity limit of up to 440 lbs. with a top speed of 4 mph. The 440 HD PLUS is available in four frame colors with five optional cushion and backrest colors that allow you to customize the color of your wheelchair fully. It's equipped with a 360-degree waterproof universal intelligent joystick for easy control. The intelligent joystick includes a power on/off button with a power indicator light, speed selection buttons, and a horn. The 12AH lithium-ion battery powers two 250W brushless DC motors that deliver a range of up to 12 miles.

26AH batteries gives up to 25 miles of range. Battery and Charger

440 HD PLUS comes in two battery options, a 288WH (24 volt, 12AH) lithium-ion or a 480WH (24 volt, 20AH) lithium-ion. 12AH battery is airline and cruise approved and delivers a range up to 12 miles. The 26AH battery is not airline or cruise approved and delivers a range up to 25 miles. The battery charging time is for each battery is 6 to 8 hours.

440 HD PLUS battery can be charged with a regular home outlet. The charging details for each battery (12AH or 26AH) are AC 100 – 240 Volt, 50 Hertz, at 1.5 – 5 amperes. 440 HD PLUS comes with a standard 24V 2A charger in the box. For quicker charging times, we recommend the 24V 3A charger for the long range 26AH battery.

The 24V 3 charger is sold separately.

Range

The highly efficient lithium-ion battery provides plenty of range.

The 12AH battery delivers up to 12 miles / 18 km. The 26AH battery delivers up to 25 miles / 40 km.

Driving habits, temperature, weather, loading weight, road environments, and other factors may affect the range.



FRAME



MOTOR & ELECTRONICS